



Presented With Compliments From
www.DirectHomeMedical.com



Built for men. Engineered for stability.

Philips Respironics GoLife for Men mask provides stability and simplicity.

PHILIPS

RESPIRONICS

sense **and** simplicity

You spoke, we listened. Introducing GoLife for Men.

Stability and simplicity. When sleep therapy providers, clinicians, and male patients told us exactly what they wanted in a nasal pillows mask, we listened...carefully. The result? The new GoLife for Men – our first sleep therapy mask built exclusively for men.

We set out to make GoLife the most stable nasal pillows mask available, starting with a no-nonsense design using facial contour arms that hug the patient's face to maintain a secure seal and stability, even when moving during sleep. GoLife for Men is designed specifically to conform to the contours of the male facial structure and features straightforward, pre-formed headgear and self-adjusting, optimally angled nasal pillows for a one-step fit.

For simplicity, GoLife has few parts, comes with all sizes, and is easy to size and fit. This can result in less time fitting the mask, and may lead to more effective titrations, fewer callbacks, and minimal to no adjustments during regular use.

Simple fitting steps

Pull the headgear over your head.




Place the cushion into your nostrils.
Select a cushion that fills your nostrils comfortably without any gaps between the cushion and the walls of your nostrils.




Adjust the side straps.
Do not overtighten.



If needed, adjust the top headgear strap. Allow the headgear to rest just above the ears like a pair of eyeglasses. The back strap will sit at the base of the head.

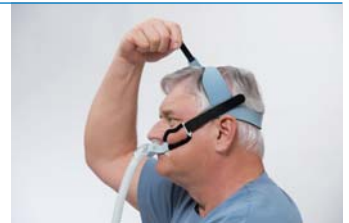

Connect the tubing to the device.
Turn on the airflow.

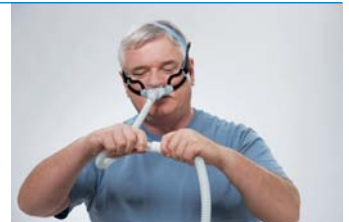


Make final adjustments while lying down.



Simple headgear is designed for fast adjustment.

Facial contour arms embrace the face with stability by hugging the cheekbones.

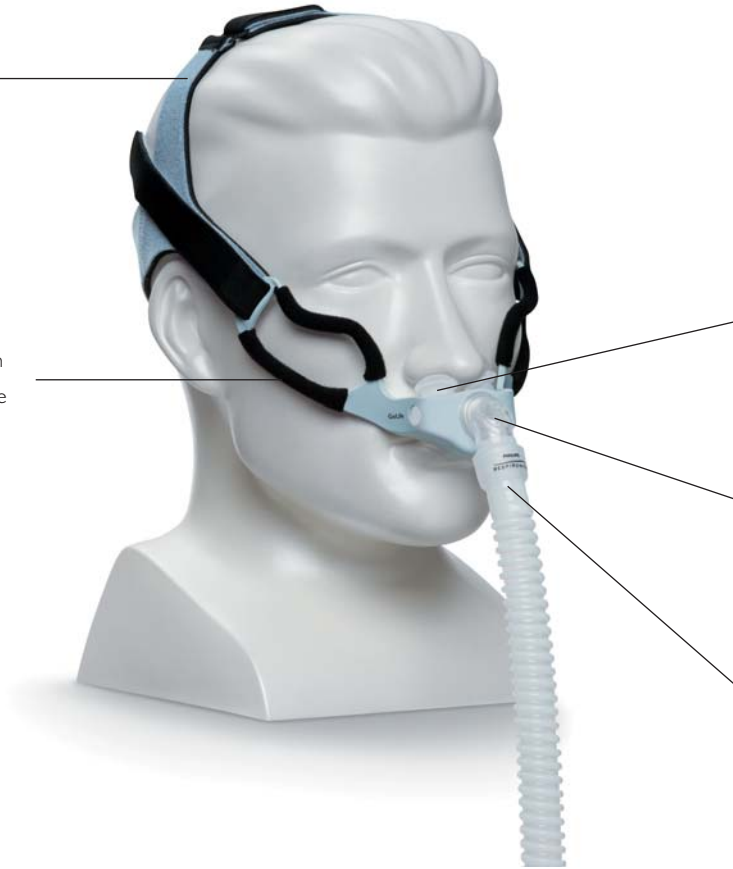

Engineered for stability. Designed for simplicity. Built for men.

Quite simply, GoLife for Men takes the complexity out of achieving successful nasal pillow therapy for male patients.

Self-adjusting, optimally angled nasal pillows eliminate the need for manual adjustment and provide a seal that's flexible.

Diffused airflow for quiet exhalation.

Long, multi-directional 18 inch/45.7 cm tube with 360 degree swivel maintains seal during patient movement.

Coming soon



Our GoLife for Women mask will offer the same customized features designed for a female patient's face.

Ordering information		
GoLife for Men mask	Part number	U.S. HCPCS code
FitPack: mask with headgear plus petite, small, medium, and large cushions	1073114	A7034, A7035
FitPack: mask only plus petite, small, medium, and large cushions	1073329	A7034
Replacement parts and accessories	Part number	U.S. HCPCS code
Mask subassembly without pillows or headgear	1073328	*
Mask subassembly (frame only, no swivel tubing)	1078408	N/A
Replacement cushion, petite	1073332	A7033
Replacement cushion, small	1073333	A7033
Replacement cushion, medium	1073334	A7033
Replacement cushion, large	1073335	A7033
GoLife for Men headgear	1073336	A7035
Tube management clip, 2 per pack	1073337	N/A
Swivel tubing (18 inch/45.7 cm) with exhalation port	1078407	A7045

*Requires addition of mask cushion for HCPCS A7034

**Philips Healthcare is part of
Royal Philips Electronics**

How to reach us

www.philips.com/healthcare
healthcare@philips.com

Asia
+49 7031 463 2254

Europe, Middle East, Africa
+49 7031 463 2254

Latin America
+55 11 2125 0744

North America
+1 425 487 7000
800 285 5585 (toll free, US only)

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)

Philips Respironics International
Headquarters
+33 1 47 28 30 82

Philips Respironics Asia Pacific
+65 6882 5282

Philips Respironics Australia
+61 (2) 9666 4444

Philips Respironics China
+86 021 24127311

Philips Respironics Deutschland
+49 8152 93 06 0

Philips Respironics France
+33 2 51 89 36 00

Philips Respironics Italy
+39 039 203 1

Philips Respironics Sweden
+46 8 120 45 900

Philips Respironics Switzerland
+41 6 27 45 17 50

Philips Respironics United Kingdom
+44 800 1300 845

www.philips.com/respironics

Respironics, FitPack, and GoLife are trademarks of Respironics, Inc. and its affiliates. All rights reserved.

Please visit www.philips.com/golife



© 2010 Koninklijke Philips Electronics N.V.
All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.
Hoech JJ 12/14/10 MCI 4103920 PN 1082047