

Oracle™ 452

ORAL MASK



Presented With Compliments From
www.DirectHomeMedical.com


Specification Sheet

An Oral Delivery Mask



The unique Oracle™ 452 delivers pressure through the mouth rather than the nose. Clinically proven to function as effectively as a nasal mask, the Oracle™ 452 is a perfect solution for mouth breathers and patients with chronic nasal obstruction.*

Oral Delivery - Oracle™ 452

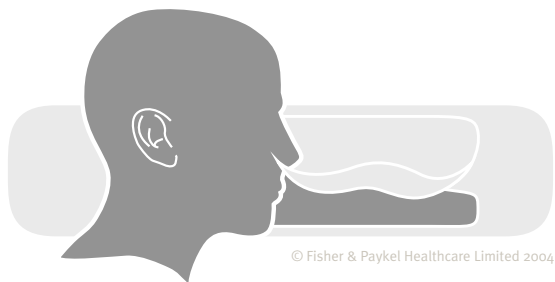

The Oracle™ 452 Oral Mask is easy to use and easy to seal. With nothing covering the field of vision, the Oracle™ 452 provides freedom of movement equivalent only to that of the Infinity™ 481 Direct Nasal Mask. The Oracle™ 452 Oral Mask comes as an Out-of-Box Solution fully assembled and ready for use.

UNIQUE MASK FEATURES FOR GREATER FREEDOM, COMFORT AND COMPLIANCE




Single Strap Headgear Minimal Intrusion

Minimal, lightweight and comfortable, single BREATH-O-PRENE™ strap allows clear vision and movement. Patient has freedom to choose Flexitube position with 360° of movement.




Rotating Adjustment Dial Easy Adjustment

Offers continuous adjustment of the distance between the SnapFlap™ Cover and SoftSeal Mouthpiece for a personalized fit. Easy to use - just rotate the adjustment dial. This can even be adjusted while mask is fitted.

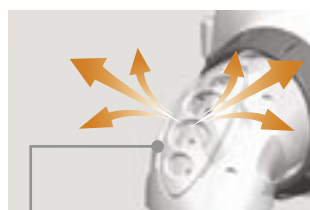

Efficient Bias Diffuser Quiet usage

Material Filter
This advanced air diffuser system prevents air from blowing onto the sleeping partner and decreases noise that could potentially disrupt sleep.

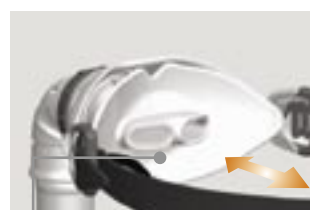

SoftSeal Mouthpiece/ SnapFlap™ Cover Interaction

Holds the lips snugly between the inner SoftSeal Mouthpiece and outer SnapFlap™ Cover providing a leak-free combination. These technologies work together to provide a simple, superior seal and maximum comfort.



Easy-Release Strap One-step removal

Simply pull the Easy-Release Cord for one-step removal of headgear. No ongoing headgear adjustments required.

CONTENTS

- Assembled mask complete with Headgear and Nasal Plugs
- Alternate size Silicone Softseal Mouthpiece included (S,L)
- User friendly Fitting Guide

REORDER CODE: HC452

REFERENCES

* Smith PL, O'Donnell CP, Allan L et al. A physiologic comparison of nasal and oral positive airway pressure. *Chest* 2003; 123 (3):689-94. * Anderson FE, Kingshott RN, Taylor DR et al. A randomized crossover efficacy trial of oral CPAP (Oracle) compared with nasal CPAP in the management of obstructive sleep apnea. *Sleep* 2003; 26 (6):721-6. * Khanna R, Kline LR. A prospective 8 week trial of nasal interfaces vs. a novel oral interface (Oracle) for treatment of obstructive sleep apnea hypopnea syndrome. *Sleep Med* 2003; 4 (4):333-8. * Beecroft J, Zanon S, Lukic D et al. Oral continuous positive airway pressure for sleep apnea: effectiveness, patient preference, and adherence. *Chest* 2003; 124 (6):2200-8.

FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER AND PAYKEL HEALTHCARE REPRESENTATIVE:

International P O Box 14 348, Panmure, Auckland 1134, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

USA Tel: 1800 446 3908 +1 949 470 3900 Fax: +1 949 470 3933 **Australia** Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 78 644 0924

Fax: +31 78 644 0914 **China** Tel: +86 20 8363 4286 Fax: +86 20 8363 3595 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7182 93777 0 Fax: +49 7182 93777 99

India Tel: +91 80 5123 6041 Fax: +91 80 5123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206

Spain Tel: +34 902 013 346 Fax: +34 902 013 379 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54 **UK (EU Authorised Representative)** Tel: +44 1628 626 136 Fax: +44 1628 626 146

Nthn Ireland Tel: 0800 132 189