

# Respironics Total™ Face Mask

**A totally unique interface solution for noninvasive ventilation and CPAP therapy**

**A Performance Series™ Mask**




The Respironics Total Face Mask is part of the Performance Series™ line which is designed to provide the right mask, the first time and every time.

The Total Face Mask is designed for patients in the hospital or home who have difficulty tolerating nasal masks using bi-level or CPAP systems. It provides several benefits:

- Unique alternative for patients who cannot tolerate nasal or full face masks
- Fast, effective seal with no fitting
- Quick release cord enabling fast mask removal
- Ideal for patients in need of immediate ventilation
- One easy-to-use size



Presented With Compliments From  
[www.DirectHomeMedical.com](http://www.DirectHomeMedical.com)

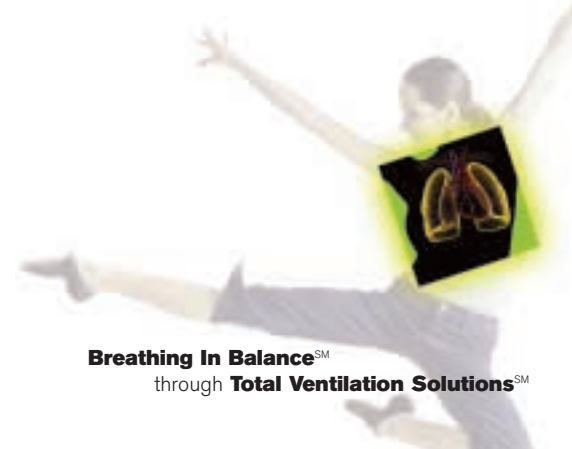


**Breathing In Balance™**  
 through **Total Ventilation Solutions™**

# Respironics Total™ Face Mask

The Total Face Mask is an interface for adult patients using bi-level or CPAP systems in the hospital or home. It is a comfortable alternative for patients who may not be able to obtain a good seal with a nasal mask, experience skin breakdown, or who are claustrophobic.

Because the Total Face Mask covers the entire face, there are no pressure points around the nose to cause sores or skin breakdown. An effective seal is created around the outside of the face as air pressure from the system inflates the soft, flexible sealing layer. This minimizes leaks while providing optimal airflow through the nasal passages. And, because the air pressure is able to circulate throughout the mask, it allows more natural, comfortable breathing for the patient. Made of clear plastic, the lightweight faceplate adds to the comfort of the Total Face Mask. It does not obstruct the patient's vision, which helps to minimize the feelings of claustrophobia that some patients experience.

The built-in exhalation device can be used in conjunction with disposable circuits to provide an additional method of exhausting CO<sub>2</sub>. The Total Face Mask comes in one size with easy-to-use headgear that can be quickly fitted on most patients – an important advantage in critical care situations. Built-in safety features, such as an entrainment valve and quick release cord, give the care provider peace of mind.

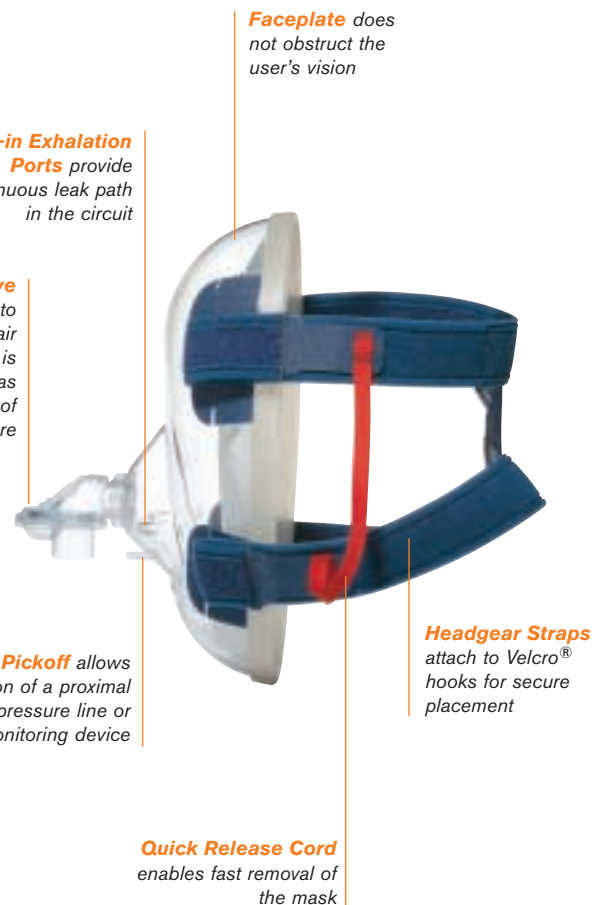

## Respironics Total Face Mask

comes with comfortable, easy-to-use headgear. Sturdy cloth sidestraps and a mesh crown strap keep the mask evenly and securely in place and help to prevent mask leaks. Adjustable Velcro tabs on the straps attach to hooks on the faceplate of the mask for stable placement.



### Ordering Information

Description	Part Number
Total Face Mask with Headgear	302433
Total Face Mask Head Strap (9 per package)	1010609
Total Face Mask Head Strap (single)	1003052



**RESPIRONICS**<sup>®</sup>  
www.respironics.com

Customer Service: 1-800-345-6443 or 724-387-4000  
Respironics EAME: +33-1-47-52-30-00  
Respironics Asia Pacific: +81-3-5280-9611

"Respironics" is a registered trademark; "Performance Series" and "Total" are trademarks; "Power Program" icons, "Breathing in Balance" and "Total Ventilation Solutions" are service marks of Respironics and its affiliates. VELCRO is a registered trademark of Velcro Industries B.V. ©2005 Respironics, Inc. and its affiliates. All rights reserved.

GEYER - RT 6/24/05 MCI 4100224

PN 1022472



SLEEP  
MANAGEMENT



CHRONIC RESPIRATORY  
MANAGEMENT



TOTAL VENTILATION  
SOLUTIONS



RESPIRATORY  
DRUG DELIVERY



NEONATAL &  
INFANT CARE