

Delivering nasal high flow with Optiflow™ at home
and in long-term care




[CLICK TO LEARN MORE](#)

Introducing the F&P myAIRVO™ 2, delivering respiratory support and airway hydration with a focus on comfort and ease of use.

The F&P myAIRVO 2 has been designed specifically to be an integrated solution for delivering nasal high flow with Optiflow. Fisher & Paykel Healthcare's leading humidification technology allows the AIRVO™ 2 to comfortably provide high flows of air/oxygen mixtures to spontaneously breathing patients, through the unique Optiflow nasal interface.


Optiflow+ Nasal Interface

- Soft, flexible prongs
- Wide-bore design reduces gas jetting
- Prongs contoured to the patient's nose

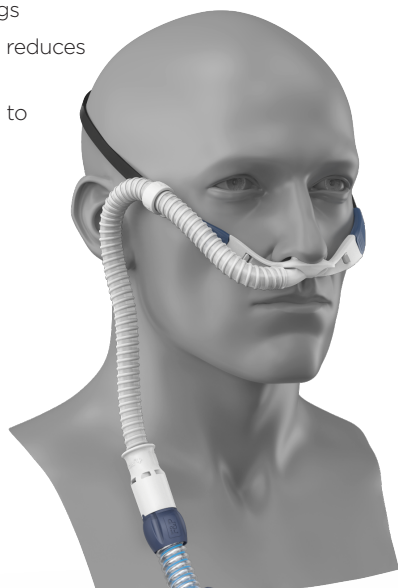

Adjustable temperature and flow settings

- Three temperature settings (37, 34, 31 °C) help achieve comfort and compliance
- Integrated flow generator delivers a wide flow range (2 - 60 L/min) - no wall air supply required

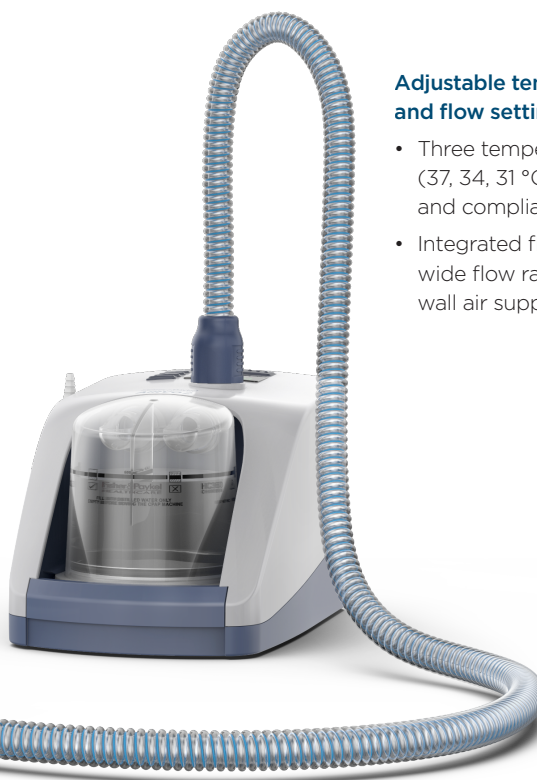

AirSpiral™ heated breathing tube

- Dual spiral heater wires and unique integrated temperature sensor
- No separate temperature probes or heater-wire adapters required

Supplementary oxygen (when required)

- Oxygen can be added from a wall supply or cylinder
- Inbuilt ultrasonic oxygen analyzer requires no calibration, service or replacement

| Product codes | |
|---------------|---|
| PT100US | myAIRVO 2 humidifier |
| OPT942E | Optiflow+ small interface (1-pack) |
| OPT944E | Optiflow+ medium interface (1-pack) |
| OPT946E | Optiflow+ large interface (1-pack) |
| OPT970E | Optiflow+ Tracheostomy interface (1-pack) |
| OPT980E | Optiflow+ mask interface adapter (1-pack) |
| 900PT560E | AirSpiral heated breathing tube (1-pack) |
| 900PT560 | AirSpiral heated breathing tube (10-pack) |
| 900PT290E | MR290 chamber (1-pack) |
| HC360 | Reusable water chamber |
| 900PT400 | Compact stand |
| 900PT401 | Water bag (2-pack) |

Optiflow, myAIRVO, AIRVO, AirSpiral are trademarks of Fisher & Paykel Healthcare Limited.
For patent information, see www.fphcare.com/ip

Find out more at www.fphcare.com/myairvo

Fisher & Paykel
HEALTHCARE