

PARI Trek® S Portable Aerosol System

	Part #	NDC #	HCPCS #
Standard	047F45-LCS	83490-0470-06	E0570, A7005
Deluxe	047F35-LCS	83490-0470-05	E0571, A7005

All configurations include Wing Tip™ Tubing



Trek® S Premium Packaging



Deluxe Carrying Case



AC Adapter & 12V DC Adapter



Trek® S Instructions & Instructional DVD



Li-Ion Rechargeable Battery
Deluxe Only



PARI LC® Sprint Reusable Nebulizer and backup reusable nebulizer





Timestrip® - 6 Month Reminder



Trek® S Compressor

The Difference is in the Delivery

Trek® S vs. Other Portable Nebulizer Systems²

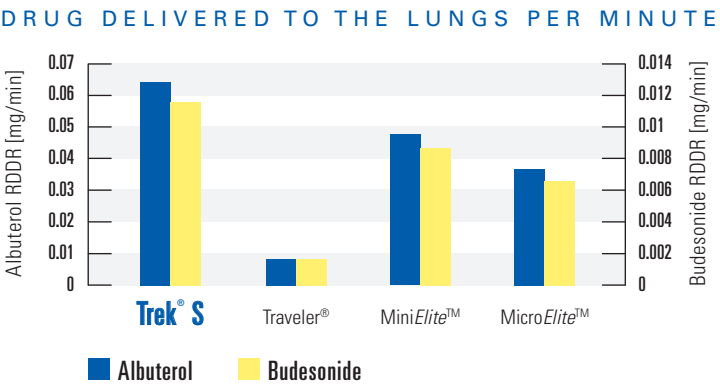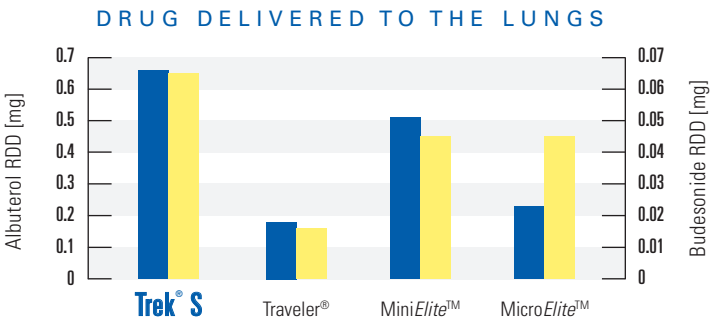

Performance with PARI LC® Sprint³

Total Output Rate (TOR) = **440** mg/min
Mass Median Diameter (MMD) = **3.8** µm
Percentage Respirable Fraction (%RF) = **65** %
Average Treatment Time (ATT) = **5-6** min

Technical Specifications:

Maximum Pressure of Compressor = **35** PSI
Maximum Liter Flow of Compressor = **8.5** LPM
Nominal Operating Pressure⁴ = **14.5** PSI (**1.0** bar)
Nominal Air Flow⁴ = **4.0** LPM
Weight = **0.8** lbs (**1.2** lbs with Battery)
Compressor Warranty = **3** years
Battery Warranty = **6** months
Battery Life = **50** min
Time to Recharge = **2** hours

Accessories:

Portable Power Pack = **047F25**
PARI LC® Sprint = **023F35**
Replacement Air Filters = **041E4851P12**



Presented With Compliments From
www.DirectHomeMedical.com



PORTABLE AEROSOL SYSTEM

*From the leader in clinically
proven aerosol delivery*



PARI nebulizers are the most widely used nebulizers by pharmaceutical companies in clinical trials for new medications compared to any other brand.



WHY TREK® S?



“Fast treatments combined with ultimate portability allow me to enjoy an active lifestyle.”



“Because patients have the freedom to take treatments anywhere, I am confident they will adhere to their treatment plan.”



“I depend on the Trek® S for safe aerosol therapy no matter where we are - at home, school, or even the big game.”

Experience the freedom of portable aerosol therapy! The PARI Trek® S Portable Aerosol System is designed from PARI's long legacy of clinically proven products and incorporates fast, efficient aerosol treatments with portable technology. Perfect for those with active lifestyles, Trek® S allows patients the freedom to take treatments anywhere.

¹Lester, et al. Nebulizer Use and Maintenance by Cystic Fibrosis Patients: A Survey Study. *Respiratory Care* 2004; 49(12): 1504-1508.

²Measured with Sierra Cascade Impactor; Albuterol 2.5mg/3mL; Budesonide .25mg/2mL; Inspiratory Flow 30liters/min; Breathing pattern (500mLs, 15cycles/min, duty cycle 0.5)

³Measured with Malvern MasterSizer X at 50% relative humidity, 0.9% NaCl solution, inspiratory flow 20 liters/minute, continuous nebulization, 23° C, fill volume 2.5 ml.

⁴When used with a PARI LC® reusable nebulizer.

^{*}All registered marks and trademarks belong to their respective companies.

Consistent Aerosol Delivery

EXCLUSIVE TO PARI, **TIMESTRIP®** IS AN EASY TO USE, VISUAL INDICATOR THAT CLEARLY INSTRUCTS PATIENTS TO REPLACE THEIR PARI REUSABLE NEBULIZERS EVERY 6 MONTHS. NO MORE GUESSWORK OR HASSLE WITH CALENDARS.



- Easy to use and activate
- Great for Patient Reminder programs
- Benefits of replacing PARI nebulizers:
 1. Consistent medication delivery
 2. Fast treatment times
 3. Reduces the risk for bacterial contamination¹