

Simplicity in Design.
Superiority in Performance.




Presented With Compliments From
www.DirectHomeMedical.com



Millennium™ M10 – Ten Liter Per Minute (LPM) Concentrator (with or without OPI*)

Respironics Millennium M10 Oxygen Concentrator delivers up to 10 LPM of oxygen reducing the delivery cost associated with 5 LPM and above oxygen patients. The M10 is engineered to reduce the cost of providing oxygen at higher liter flows. The M10 has fewer parts than any other concentrator sold today. Fewer parts, combined with greater reliability, help you to reduce expenditures on accessories and unnecessary inventory.

Millennium M10: Design and Performance

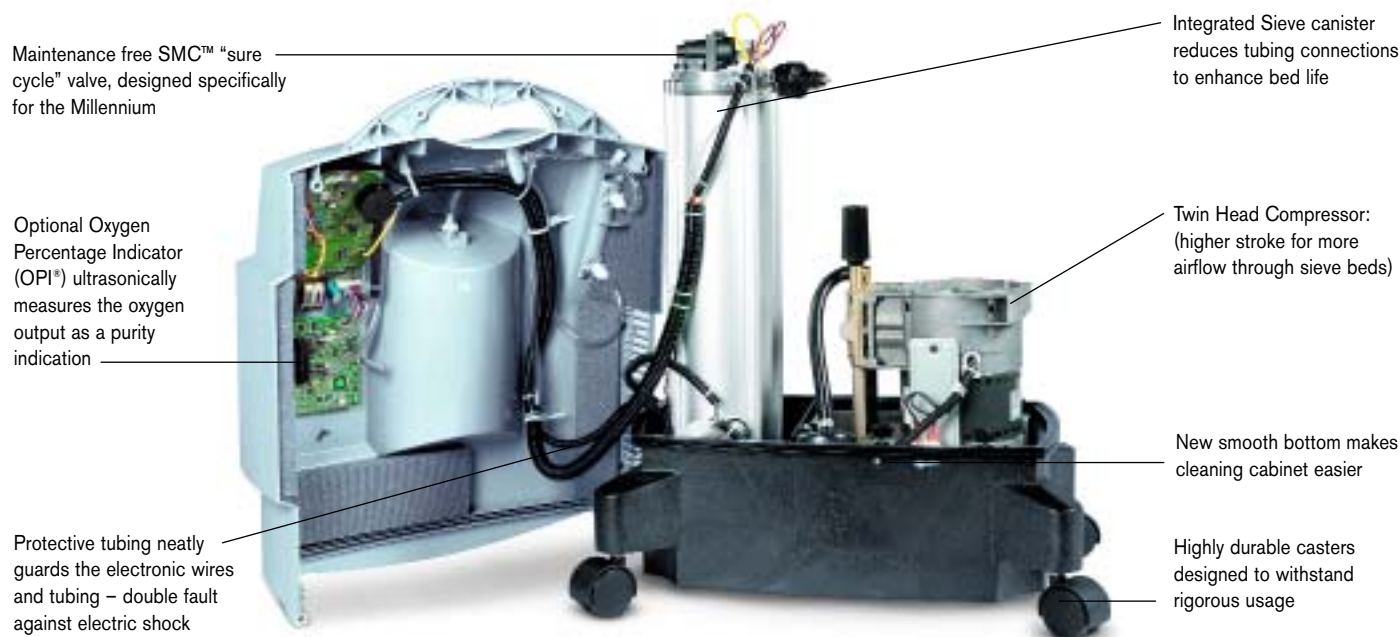

Specifications

Flow Rate	10 LPM in 1 liter increments
Electrical Requirements	120V/60Hz (± 10%)
Oxygen Percentage Indicator	Green Light greater than 82% Yellow Light between 70-82% Red Light less than 70%
Oxygen Concentration*	92 +/- 4% @ 8-10 LPM 94 +/- 2% @ 3-7 LPM 92 +/- 4% @ 1-2 LPM
Weight	53 +/- 0.5 lbs
Dimensions	27" x 19" x 13"
Storage/transport temperature	-30 to 160° F
Operating temperature	55 to 90° F
Storage/transport humidity	Up to 95%, non-condensing
Operating pressure	10-30 psig
Warranty	One year standard Optional two years (excludes compressor) Optional two year compressor
Alarm indicators	High system pressure Low system pressure Power failure Low oxygen level (OPI® model available on the M10605 model only) No oxygen flow

* Device operation above or outside of the voltage, LPM, Temperature, Humidity and/or altitude values specified may decrease oxygen concentration levels.

Covered by one or more of the following U.S. Patents: 5,060,506; 5,183,483; 5,916,349; 5,997,617; 6,190,441; 6,348,082; 6,382,931; 6,395,065; and 6,497,755

Respironics' Millennium M10 Oxygen Concentrator is designed to be the highest performing and most reliable machine available today. Each M10 concentrator is individually tested in accordance with Respironics' quality guidelines and comes with a written certificate of conformance and one year standard warranty.

The optional Oxygen Percentage Indicator® (OPI) is an audible alarm that sounds when oxygen levels fall below 70%.

Ordering Information

Description	Part Number
Millennium M10 without OPI	M10600
Millennium M10 with OPI	M10605
Accessories	Part Number
25' Tubing, Connector, and High Flow Cannula	1014653



Customer Service: 1-800-345-6443 or 724-387-4000
 Respironics Europe: +33-(0)1-55-60-19-80
 Respironics International: +33-(0)1-55-60-19-80
 Respironics France: +33-(0)2-51-89-60-00
 Respironics Deutschland: +49-81-52-930-60
 Respironics Asia Pacific: +852-234-342-18

CAUTION: U.S. federal law restricts this device to sale by or on the order of a physician. Millennium and OPI are trademarks and "People.Products.Programs" is a servicemark of Respironics, Inc and its affiliates.

©2004 Respironics, Inc. and its affiliates.

Gyer KB 3/29/05 MCI 4100221 PN 1023229



SLEEP
MANAGEMENT



CHRONIC RESPIRATORY
MANAGEMENT



TOTAL VENTILATION
SOLUTIONS



RESPIRATORY
DRUG DELIVERY



NEONATAL &
INFANT CARE